


# INSTRUCTIONS FOR Warming Cabinets (#4833, 4834, 4857 & 4858 )



## Key Features

- Touch-safe heater with a protective shield to minimize burn risk.
- All components are UL-recognized for safety and reliability.

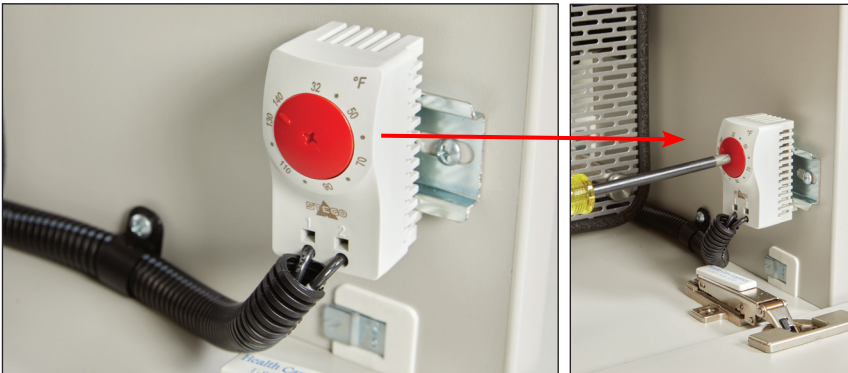

**Step 1.** Use a Phillips screwdriver to preset cabinet to the desired temperature. Close the door and allow the unit to warm for 24 hours. The adjustable thermostat allows control from room temperature up to 140°F ( $\pm 7^\circ\text{F}$ ), while the quiet, low-speed fan evenly distributes heat.



**Step 2.** After 24 hours place content inside of the warming cabinet and shut the door.

**NOTE:** The end user is responsible for setting the thermostat and monitoring temperatures with a data logger to ensure there are no excursions outside the drug manufacturer's specified requirements.